


HOT & COOL
SKIN FST

Product Configuration

Product Configuration



Hot & Cool Skin Fit body

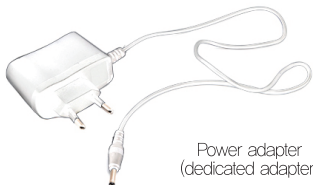

Power adapter
(dedicated adapter)



Cradle



User manual

This product offers hot and cold massage features which assist in maintaining healthy skin.
Please read the user manual carefully for proper and safe use.

Product Description

Product Description



Model	HCS-1000	Dimensions(mm)	45,5(W)x51,5(L)x140(H)
Weight	126g	Color	White / Orange
Power input	Adaptor : 100V~240V (5V 500mA)	Battery	Li-Polymer 1300mA Rechargeable
Temperature	Heating : 39°C / Booster : 36°C / Cooling : 13°C / Tightening : 5°C		


* In Room Temperature Standard, Error $\pm 2^{\circ}\text{C}$, It can be different by user's environment

How to use

Basics

- Remove the protective tape from the head before using for the first time.
- Press the power button to turn on the device and enter the heating mode.
- You can select between four different modes by pressing the Power/Mode button.
- Each mode will run for 3 minutes.
- The device will enter the standby mode when each mode is completed.
The LED for the next mode will blink while the device is in standby mode.
- Press the Power/Mode button in standby mode to enter a specific mode.
- The device will automatically switch off when left for 5 minutes.
- Press the Power/Mode button for more than 1 second to turn off the device.
- Always clean the head and device with a soft cloth, cotton pad or swab after use.

How to use the Power/Mode button

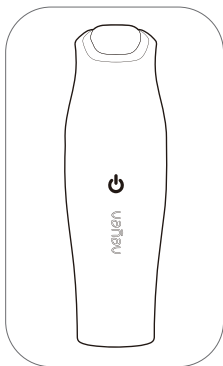

- Press the Power/Mode button as described below to enter a specific mode.

1 time	Power ON / Heating mode / Red LED
2 times	Booster mode / Yellow LED
3 times	Cooling mode / Purple LED
4 times	Tightening mode / Blue LED

- Pressing the Power/Mode button 1~4 times will change the mode as described above.

- * Press the Power/Mode button for more than 1 second to turn off the device.
- * It can feel cosmetic absorption effect deeply with micro vibration in
Booster / Tightening Mode

How to charge

How to charge Hot & Cool Skin Fit

- This device can be used as a portable device after charging.
- * The device cannot be used while charging. Please disconnect the power adapter before use.



Connect the power adapter to the connector located on the bottom of the device.



Connect the power adapter to a power outlet.



Put the device in the cradle as described in the image and run the cord out from behind.



Charging
Red
LED



Charging complete
Green
LED

How to use each mode

Facial care



Enter Heating mode and gently massage around the skin.



Apply moisturizer or nourishing cream after initial basic skin care.



Enter Booster mode and massage the skin along the grain.



Finish with Tightening mode after massaging.

When using a mask pack



Enter Heating mode and gently massage around the skin.



Apply the mask pack.



Enter Booster mode and massage the face with the mask pack on.



Finish with Tightening mode after massaging.

When using a nose pack



Enter Heating mode and gently massage around the nose.



Apply the nose pack.



Enter Tightening mode and massage after removing the nose pack.

When applying makeup



Enter Tightening mode and massage the skin pores before applying makeup.



This will temporarily reduce the pores and offer porcelain-like skin after makeup.

Caution


Caution before use

- Do not use on children or impaired persons.
- Please consult a doctor before using during pregnancy.
- Do not use if there is a skin problem including wound, eczema or dermatitis.

Safety precautions

- Please read carefully before use.
- The safety precautions described in this section include important contents related to safety. Please follow the directions carefully.
- This manual categorizes risks associated with misuse and wrong handling into three categories.

 **Danger!** Negligence of this symbol can lead to death or severe injuries.

 **Warning!** Negligence of this symbol can lead to severe injuries.

 **Caution!** Negligence of this symbol can lead to injuries or damage.

 **Danger!** *This device should not be used on people who are using one or more of the following devices.*

- Implanted electronic devices including pacemakers.
- Life-supporting electronic medical devices including heart-lung machine.
- Attachable medical electronic devices including ECG monitor.

 **Warning!** *This device should not be used on people with the following symptoms or on the following body parts.*

- Metal allergy / injury / impaired pain sensation or perceptual disturbance / severe heart disorder / fever / children under 12 years of age
- Severe blood pressure abnormality / pregnancy / bleeding / dermatitis / infectious diseases / ear / mouth (uvula) and genitals
- Fractured part / irritated skin due to excessive tanning / inflammation caused by acne

**Caution!**

*Do not drop or apply shock to the device. This can cause damage.
The user will be charged for repairs for defects caused by the user even during the warranty period.*

- Turn off the device and connect to the power adapter when not in use.
- Disconnect power when the device is not used for seven days or more, and check safety before using again.
- Charge the device regularly before the battery is completely depleted. Frequent depletion of the battery can cause deterioration of the device's performance.
- If you have sensitive skin, test the product on the inside of the forearm to check for adverse reactions.
- Do not press too hard when massaging. This can cause skin irritation.
- Stop using the product if you feel that the temperature is too hot or cold, and please contact our customer service.
- Wash the face with warm water before using this device.
- Do not use this device for purposes other than those described in this manual.
- Do not use or store in places with high temperature, dust, humidity or fire such as the bathroom, kitchen or inside a car which can shorten the product's lifespan and cause defects.
- Use only the power adapter included in this product.
- Do not place heavy objects, water or beverage on the device.
Liquids inside the device can cause malfunction.
- Do not clean the device with ignitable materials such as alcohol, benzene or paint thinner.
- Do not disassemble or self-repair device.
- Keep out of reach of children.
- Do not use the device on a single spot for a prolonged time.
- Do not touch the power cord with wet hands when pulling the plug to disconnect the power adapter.
- Keep the device dry. If the device is dropped in water, take it out immediately and dry the device sufficiently and have it inspected at our customer service center.
- Wipe off residual cosmetics and clean the device with a soft cloth, cotton ball or swab after use.

FAQ

Q. The device doesn't work.

- A The device may not work if the battery is depleted.
Connect the device to the power adapter for charging.

Q. The power button doesn't work when charging.

- A The device will not work during charging.
Please use the device after disconnecting it from the power adapter.

Q. The device beeps after the blue and red LED blinks, and the power turns off automatically.

- A This can occur if the cooling fan isn't working when the battery is depleted.
Charge the device before use. If the symptom continues to occur after charging, the cooling fan may be damaged or is malfunctioning due to a foreign object.
Please contact our customer service center.

Q. The head becomes hot again after using the cooling mode /tightening mode.

- A The heat remaining after using the cooling mode is being dissipated through the head.
This is normal and the head will cool down over time.
Please use the device again after a couple of minutes.

Q. The Power/Mode LED is blinking.

- A The device offers 4 modes and each mode runs for 3 minutes.
The LED for the next mode will blink after each mode is completed after which the device will enter Standby mode. The device will automatically switch off when left for 5 minutes.

Q. The charging LED is not blinking and the device is not charging.

- A Check if the device is properly connected to the power adapter.
Charging may not occur if the power adapter's connector is dirty.
Clean the power adapter's connector and connect the charger again.

Q. The power went off on its own during use.

- A The power may go off on its own if the battery is low.
Connect the device to the power adapter and charge sufficiently before use.

Q. Device doesn't work Micro vibration.

- A Micro vibration works only in Booster / Tightening mode. If Micro vibration doesn't work in Booster / Tightening mode, please contact our customer service center.

Q. I have a problem that is not described above.

- A Please contact our customer service center.

