

• Safety Precautions

It is not intended for use by children, pregnant women or patients with heart disease. Please do not use it on wounds, reddish and swelling skin and skin with eczema.

• Danger

Do not use this device if you are using any of the following medical devices:

1. Implanted medical device e.g. cardiac pacemaker
2. Artificial organ e.g. heart-lung machine
3. Implantable electronic medical device e.g. electrocardiogram

• Warning

Do not use this device on the following patients or areas: pregnant women, patients with allergy to metal, severe heart disease, fever, abnormal blood pressure, infectious disease, or areas with trauma, hemorrhage, dermatitis, mouth (uvula) or vagina, or fracture.

• Caution

1. Do not use the massager for any purposes other than intended as it may cause accident.
2. Do not place any heavy object or water (liquid) near the device. Any water or liquid inside the device may cause damage.
3. Do not disassemble or repair the massager without authorization.
4. Do not drop or apply impact to the massager.
5. Make sure to wash your face and remove any dirt before using the device.

※ If the LED lamp does not turn on or the vibration is too weak, make sure to replace the batteries.

※ If the massager does not operate properly when making contact with skin, make sure to replace the batteries.

Directions for Use

1. Place AAA alkaline batteries into the main unit.
2. Apply a cosmetic product to your face.
3. While holding the touch sensing electrode panels located on each side of the main unit, press the head part of the massager against the area you want to massage.
4. Gently massage your face along the skin texture with the head part of the main unit.
5. Keep it clean using tissue, soft cloth, pad or cotton swab after use.

※ Apply a sufficient amount of cosmetic product, otherwise the massager may not slide on skin smoothly.

※ Make sure that your hand or skin is in contact with the massager when using it on other persons.

Eye Care


1. Apply an eye care product around your eye areas.
2. Press the head against the eye areas and gently massage it as shown in the image.
3. Start from under eye areas and circle out the rest. Repeat this 2-3 times.

Recommended Skincare Products: eye cream, eye gel and eye care products


Normal Skin Care

1. Apply skincare products after washing your face.
2. Press the head against your face and gently massage it as shown in the image.
3. Start from the center of your face and massage in a spiral motion. Repeat 2-3 times.

Recommended Skincare Products: skincare products according to your skin type

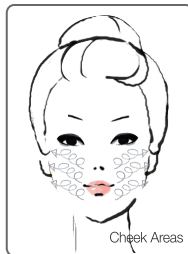


How to Replace Batteries

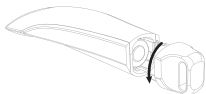

- ① Rotate the battery cover counter-clockwise to remove it.



- ② Place batteries with the positive side facing towards the front part of the massager.

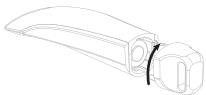

- ③ Rotate the battery cover clockwise to close it.

Caution

- ※ Do not use excessive force to rotate the battery cover in the wrong direction.
- ※ Make sure to use AAA alkaline batteries.

FAQ

Q. The massager and its LED lamp do not turn on.

A. Make sure to hold the touch sensing panels located on each side of the main unit or check the batteries and replace them if necessary.

Q. Can I only use it on my face?

A. The massager is designed to massage your face. However, it can be also used on other parts of your body including hands and neck.

Q. The massager does not slide on my face when a cosmetic product is fully absorbed. What should I do?


A. The massager may not slide smoothly on skin if you do not apply a sufficient amount of cosmetic product. In this case, you may apply more product depending on your skin condition or tap it lightly on your skin.



Company : Doctors Tech Co., Ltd.
Product Type : Cosmetic Device (Skin Massager)
Model : Ray-1001
Manufacturer : Doctors Tech Co., Ltd.
Country of Origin : Made in Korea

uānau

RAY-1001

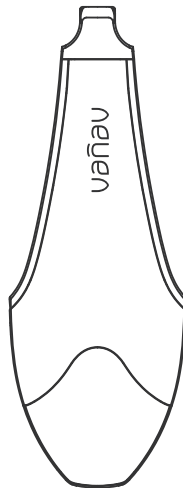

Please read the user instructions carefully before use.